


Backing up to a USB Memory Stick

First, go to the System Tab and then refer to the image below.

1. Make sure the checkbox (Automatically backup my data at the end of each session) is checked
2. Click the button (...) at the right edge of the Backup Location text box

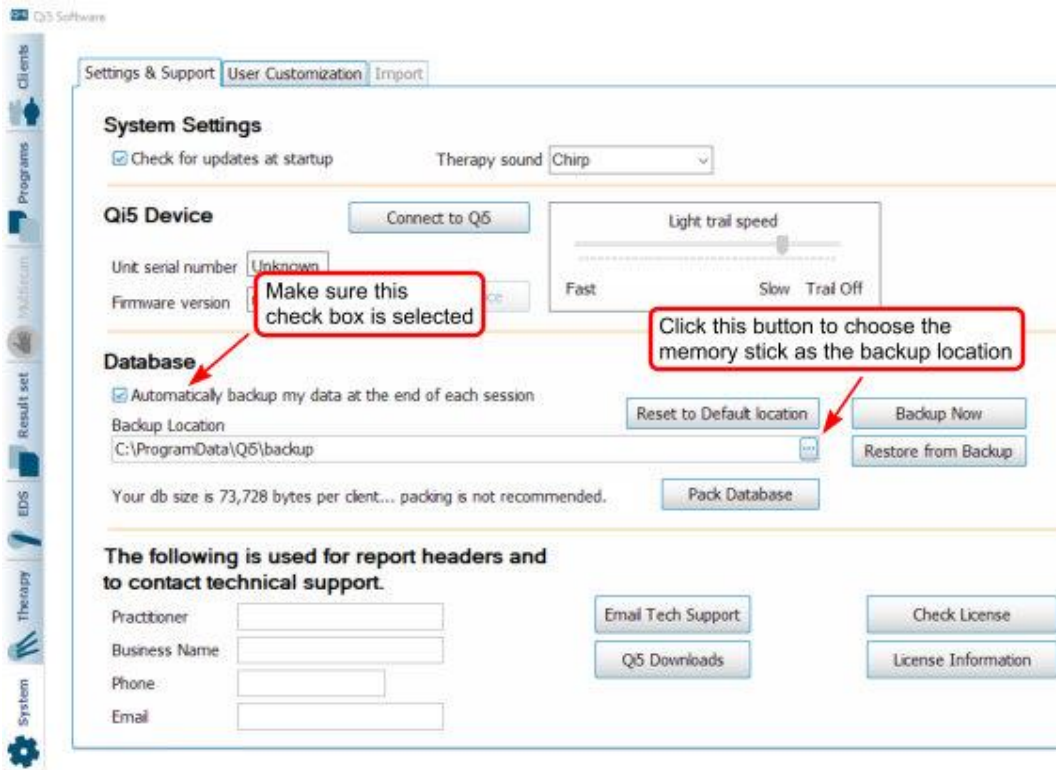

When you click the (...) button a popup will appear which allows you to select the USB memory stick location. You need to select the USB memory stick drive using this popup. The name and letter of the USB memory stick can be different on every system. On my system, the popup looks like this...



Click on the drive representing the USB memory stick and click the OK button. When you return to the System screen, you will see that drive name showing in the Backup Location text box.

Now just close the Qi5 software and your database will be compressed and written to the USB memory stick in a zip file. Note that the backup location will remain until you change it again, so if you remove the USB memory stick, the software will no longer be able to make a backup. You can use the button labeled (Reset to Default) on the system screen to set the backup location back to its original setting.